


GEOSTATS PTY LTD

Mining Industry Consultants
Reference Material Manufacture and Sales

Assay Pill Product Code

GAP-02

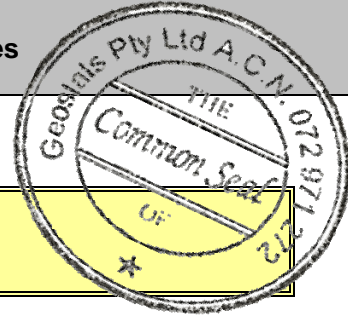

Gold grade of 1g pill (ppm):	1025
Gold grade of one pill in 1kg of blank (ppm):	1.025

Assay Pill Details

Preparation and Packaging

These assay pills are made by precipitating gold from a weak aqua regia solution. Appropriate amounts of silica and cement are blended to attain the target grade. 1g of the powder is pressed to form each pill. 25 x 1g pills are packed into glass jars ready for purchase and use.

Usage

Pills should be inserted into a coarse blank material with typical physical characteristics. The mass of the blank material needs to be recorded so results can be assessed. Multiple pills can be used if higher grades are required. Conversely, more blank material will lower the expected grade. Visit our website for a useful Excel spreadsheet for determining grades.

Assay Testwork

Two pills are analysed by LabWest in Malaga to confirm the target grade has been achieved. The value attributed to the pills is by mass calculation.

10A Marsh Close, O'Connor, Western Australia 6163

Phone : +61 8 9314 2566, Fax : +61 8 9314 3699

e-mail : pjh@geostats.com.au, srr@geostats.com.au

Website <http://www.geostats.com.au>

GAP-02

Geostats Pty Ltd, Assay Pills Product Code :